

SONY



UP-DR80MD

Digital Color Printer



Sony is proud to introduce a new stylish digital color printer, the UP-DR80MD.

Demand for color printing in the medical environment is ever increasing, especially for A4/letter-sized prints, as more and more color images are required in medical reporting and records. The expectation is also for durability and high quality, as these prints need to be kept in long-term storage.

In response to these demands, the UP-DR80MD incorporates Sony's dye sublimation technology, which has already been well received for its high-quality prints in the digital photo market. Thanks to this technology, the UP-DR80MD produces high-quality and highly durable prints in A4/letter size with a fast printing speed.


The UP-DR80MD boasts a compact and sophisticated design with a stylish white color that blends well with medical modalities such as cart systems.

Moreover, with a wireless printing capability provided by the optional UPA-WU10 Wireless Print System, you can achieve cable-free connection between a medical imaging modality and an 'off-cart' printer without requiring a lengthy USB cable.

Thanks to its superb picture quality and high durability, the UP-DR80MD is ideal for creating full-color medical images on A4/letter-sized prints.



For open surgery

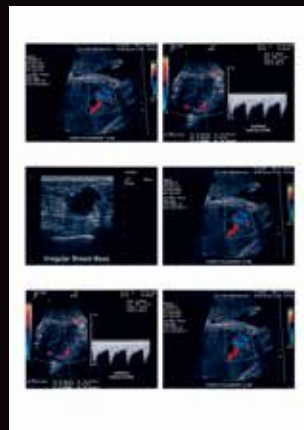

For ultrasound



For dental

Features

Superb Print Quality

Incorporating Sony's dye sublimation technology, the UP-DR80MD produces photo-quality prints. It also employs Sony's laminating technology, which coats each print with a protective layer that's resistant to light, heat, water, humidity, and fingerprints. This laminating technology allows the prints to be placed in long-term storage without fear of degradation. The UP-DR80MD can print onto two different sizes of paper: A4 (297 x 210 mm, 11 3/4 x 8 3/8 inches) and Letter (279 x 216 mm, 11 x 8.5 inches).

Stylish Design and Compact Size


The UP-DR80MD has a stylish design that fits perfectly with most medical environments. Its glossy white surface allows users to wipe it clean easily. Also, the compact body of the UP-DR80MD (approx. 317 (W) x 207 (H) x 425 (D) mm (12 1/2 (W) x 8 1/8 (H) x 16 3/4 (D) inches) enables it to fit into both medical carts and constrained spaces.

Wireless Printing Capability*

By using the optional UPA-WU10 Wireless Print System, you can achieve cable-free connection between the UP-DR80MD and surgical endoscopes, ultrasound systems, and other hospital imaging modalities with compatible wireless capability.

This is ideal for busy clinical environments, where modalities are constantly being moved between different OR and treatment rooms. This all-in-one wireless printing solution also resolves the problem of limited space on medical carts, allowing connection between a medical imaging modality and an 'off-cart' printer without requiring a lengthy USB cable.

Use Case (Operation Room)

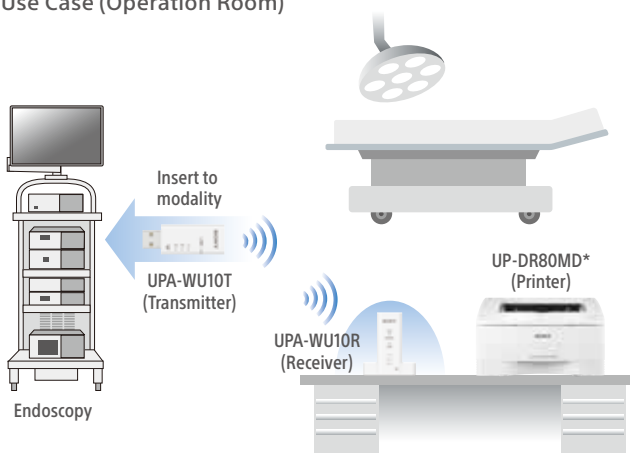

* Function available only on models with USB 5V Output port, manufactured in 2016 or later. Regarding the compatibility for wireless use, please contact your nearest Sony office or an authorized dealer.

Useful Front Operation

The front-loading operation makes it easy to replace paper and ribbon, and allows the UP-DR80MD to be installed in even the tightest of spaces. This useful design feature also allows you to use the printer on-cart.



High-speed Printing

The UP-DR80MD operates at high speed. Printing an A4-sized image takes approximately 76 seconds and a letter-sized image takes approximately 72 seconds.

Easy Color Adjustment

The UP-DR80MD offers easy color adjustment via its printer driver property. This incorporates a preview image from which you can confirm and adjust images easily.



Color Look Up Table

Two types of color look up table are available to reproduce suitable colors for different modalities, such as endoscopy and ultrasounds.

The default "Natural" setting creates a color based on the s-RGB monitor, which is suitable for endoscopy images. For ultrasound images, a "CG" color mode is available that reproduces more vivid colors. These settings are selectable via the printer property.



Gray Balance Adjustment

A Gray Balance Adjustment function allows you to easily select the color that most closely matches your requirements and provides the most suitable tonal range for your application.



Specifications


UP-DR80MD	
Printing Method	Dye Sublimation Thermal Printing
Resolution	301 dpi
Gradations	8 bits (256 levels) processing each for yellow, magenta and cyan
Effective Print Pixels	A4 size UPC-R80MD: 3400 x 2392 dots Letter size UPC-R81MD: 3192 x 2464 dots
Print Size	A4 size: 287 x 202 mm (11 1/4 x 8 inches) Letter size: 269 x 208 mm (10 5/8 x 8 1/4 inches)
Printing Time	A4 size: Approx. 76 sec Letter size: Approx. 72 sec
Interface	Hi-Speed USB (USB 2.0)
Paper Tray Capacity	50 sheets
Power Requirements	AC 100 to 240 V, 50/60 Hz
Input Current	3.4 to 1.4 A
Output Current	0.5 A/5 V (for UPA-WU10*)
Operating Temperature	5 °C to 35 °C
Operating Humidity	20% to 80% (no condensation allowed)
Storage and Transport Temperature	-20 °C to 60 °C
Storage and Transport Humidity	20% to 80% (no condensation allowed)
Dimensions*1	Approx. 317 (W) x 207 (H) x 425 (D) mm (12 1/2 (W) x 8 1/8 (H) x 16 3/4 (D) inches)
Mass	Approx. 11.5 kg (25.3 lbs)
Supplied Accessories	Paper Holder (pink) (1) Paper Holder (blue) (1) Cleaning Ribbon (1) USB Cable 1-824-211-41 (SONY) (1) Before Using this Printer (1) Quick Guide (1) CD-ROM (1) Software License Agreement (1) Service Contact List (1) European Representative (1)

*1 The values for dimension are approximate.

This product is distributed to the US and EU as a medical device and satisfies product safety standards (e.g. IEC 60601-1).**

** For more details, please contact your nearest Sony sales office or an authorized dealer.

Dimensions

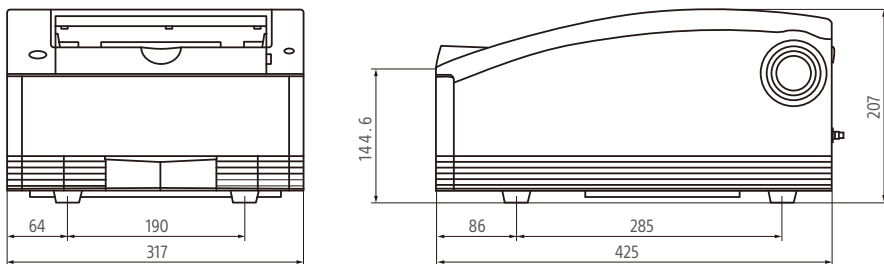

Unit: mm

Print Media

UPC-R80MD: A4 size

Contents: Paper roll for 50 prints x 2
Ink ribbon x 2



UPC-R81MD: Letter size

Contents: Paper roll for 50 prints x 2
Ink ribbon x 2



Wireless Print System (Optional accessories)

UPA-WU10*

Consists of UPA-WU10T (Transmitter) and UPA-WU10R (Receiver)



* To check region-compatibility for this product, please contact your nearest Sony office or an authorized dealer.

Distributed by

©2017 Sony Corporation. All rights reserved.

Reproduction in whole or in part without written permission is prohibited.

Features and specifications are subject to change without notice.

All non-metric weights and measurements are approximate.

Sony is a registered trademark of Sony Corporation.

All other trademarks are the property of their respective owners.

Please visit Sony's professional website or contact your Sony representative for specific models available in your region.